

Rapid SARS-CoV-2 Antigen Test

- Summary

TEST KIT
EXPRESS



Box Size Options:

- Boxes of 20 tests
- Boxes of 5 tests

- ⊕ Outperforms other Covid-19 Test Kits for **Accuracy AND Speed**
- ⊕ **Near 100% sensitivity** on highly infectious people and **99% on moderately infectious.**
- ⊕ **100% Specificity** in over 10,000 samples meant no misdiagnosed 'false positives'
- ⊕ Visual results observed within **15 minutes**
- ⊕ **CE & UK Compliant** and **MHRA registered**
- ⊕ **Developed in the UK** by **Professor Sir Chris Evans, OBE** and his Excalibur Health Team
- ⊕ **Ethically Priced** – from only **£6.25 per test (volume discounts available)**
- ⊕ Test Kit Express donate £1 for every 20 tests to **Wooden Spoon the children's charity of rugby**

*Test currently MHRA registered for professional use only

*Wooden Spoon Registration No 326691 (England & Wales) and SC039247 (Scotland)

Ethically Priced To Be As Affordable As Possible

Test Kit Performance Summary

- (Independently Tested)

TEST KIT
EXPRESS



Outperforms Other Tests

- For both **ACCURACY** & **SPEED**



Most Reliable

- Near **100% 'Specificity'** means **no misdiagnosed 'false positives'**



Most Accurate

- **100% 'Sensitivity'** for highly infectious individuals and 99% for moderately infectious individuals means **detection in almost every single case**



Fastest

- **Results in 15 Minutes.** A positive result may be seen more quickly and will be valid



Straightforward

- **Easy to Administer** nasal swab means **individuals won't mind being tested regularly**



Cost-Effective & Ethically Priced

- **From just £6.25 per test**

Test Kit Express / Excalibur Test Independently Tested by:



Ethically Priced To Be As Affordable As Possible

Background to Test Kit Express Pricing

- Why the Test Kits are so affordable

TEST KIT
EXPRESS

- ✓ Organisations and individuals can take great comfort knowing the Test Kits were developed by **Professor Sir Chris Evans, OBE** and his **UK-based** Excalibur team.
- ✓ Sir Chris is **Europe's foremost leader** in the fields of medicinal bioscience, molecular microbial testing, clinical diagnostics and pharmaceutical development.
- ✓ Having already supplied the UK and other governments with over 50 million tests, new Rapid Test Kit **outperforms other tests for both accuracy & speed.**
- ✓ Sir Chris has kept the **price of the Test Kits very low** to make them affordable to businesses, the NHS, care services and schools.
- ✓ **Test Kit Express** are part of a leading UK-based supplier of PPE, **PPE Boxes Direct** www.ppeboxesdirect.co.uk
- ✓ Following Sir Chris' ethical lead, Test Kit Express donate £1 for every 20 tests to charity: **Wooden Spoon – the children's charity of rugby (Bristol Bath Somerset region)**



Professor Sir Chris Evans, OBE

Ethically Priced To Be As Affordable As Possible

A bit more info about Professor Sir Chris Evans, OBE



Professor Sir Chris Evans, OBE

- ❑ **Founder and Chairman** of Excalibur Healthcare Services Ltd.
- ❑ Responsible for developing and manufacturing **millions of clinical diagnostic tests for global distribution.**
- ❑ Invented series of rapid tests for microbial detection. **Over 1.5 billion units of these rapid bioluminescent tests distributed.**
- ❑ Sir Chris and the Excalibur testing team continue to be at the very forefront of the industry having developed this highly reliable and accurate SARS-CoV-2 Rapid Antigen test that **outperforms other marketed Antigen tests.**

Ethically Priced To Be As Affordable As Possible

More Reasons To Be Reassured

TEST KIT
EXPRESS

- Compliance, Accreditation & Independently Tested

- ✓ **Accreditation** from multiple institutions, e.g. **BSI** (British Standards Institute), **ISO** (International Organisation for Standardization) and **TUV** (Germany's Technical Inspection Association).
- ✓ **Extensive testing** by **Imperial College London**, **Liverpool School of Tropical Medicine (WHO-backed laboratory)** and **AstraZeneca** plus many international tier 1 laboratories showed Sir Chris' Test Kit **outperforms other tests for both Accuracy & Speed.**
- ✓ The tests are **currently being used across Europe** in Germany, Italy, Austria, Belgium and the UK in **hospitals, care homes, corporate settings, sporting and educational establishments.**



Ethically Priced To Be As Affordable As Possible

Ordering & Prices

– Boxes of 20 Tests and Boxes of 5 Tests

- ✓ Ordering **online** is quick & easy



- ✓ Enquire by **phone** for more than 100 boxes



Boxes of 20 CoV-2 Rapid Antigen tests

No. of Boxes	Price Per BOX	Price Per TEST	TOTAL Delivery Charge (Secure Delivery)
1-100	£125	£6.25	£10 to £50
101+	Contact Us		

Boxes of 5 CoV-2 Rapid Antigen tests

No. of Boxes	Price Per BOX	Price Per TEST	TOTAL Delivery Charge (Secure Delivery)
1-200	£40.00	£8.00	£7 to £50
201+	Contact Us		

Ethically Priced To Be As Affordable As Possible

Box Contents

– Boxes of 20 Tests and Boxes of 5 Tests

Test Kit Box Contents – Boxes of 20 Tests

ITEM	QUANTITY
SARS-CoV-2 Antigen Test Cartridge	20
Extraction Solution Bottle	2
Extraction Tube	20
Nasal Swab	20
Instructions for Use	1

Test Kit Box Contents – Boxes of 5 Tests

ITEM	QUANTITY
SARS-CoV-2 Antigen Test Cartridge	5
Extraction Solution Bottle	1
Extraction Tube	5
Nasal Swab	5
Instructions for Use	1



Information	Details
Time to Result	Within 15 Minutes
Storage	4-30 Degrees C
Shelf Life	18 Months
Specimen Type	Nasal Swab

Ethically Priced To Be As Affordable As Possible

Specimen Collection

Training in specimen collection is recommended because of the importance of specimen quality. Use the swabs provided in the test Kit.




1. Carefully insert the swab into the nostril of the patient, reaching the surface of the posterior nasopharynx that presents the most secretion.
2. Swab over the surface of the posterior nasopharynx. Rotate the swab several times.
3. Withdraw the swab from the nasal cavity.

Ethically Priced To Be As Affordable As Possible

Specimen Preparation

1. Add 8 drops (about 0.25mL) of extraction buffer into the extraction tube.
2. Place the swab with the specimen into the extraction tube. Roll the swab three to five (3-5) times. Leave the swab in the extraction buffer for 1 minute.
3. Pinch the extraction tube with fingers and remove the solution from the swab as much as possible. Dispose of the used swab in accordance with your biohazard waste disposal protocol.
4. Install the nozzle cap onto the sample extraction tube tightly. Use the extraction solution as test specimen.

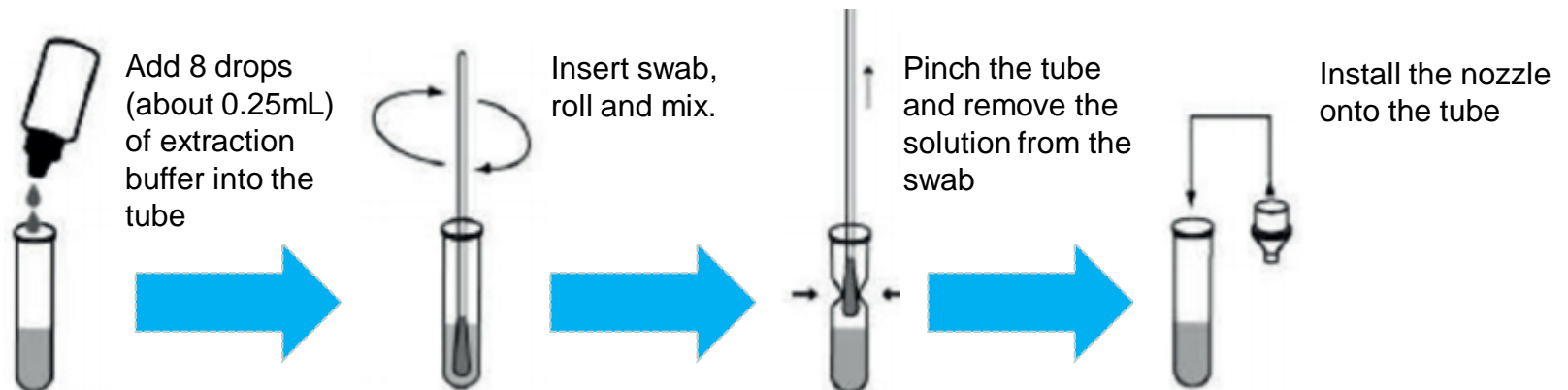


Ethically Priced To Be As Affordable As Possible

Results Procedure & Interpretation

1. Open the pouch and remove the card. Once opened, the test card must be used immediately. Label the test card with the patient identity.
2. Invert the extraction tube and add 2-3 drops (50-75 μ L) of test specimen to the well (S) by gently squeezing the extraction tube.
3. Read the results in 2-15 minutes – a positive result may be seen quickly and will be valid

Note: Results after 20 minutes should be considered invalid.

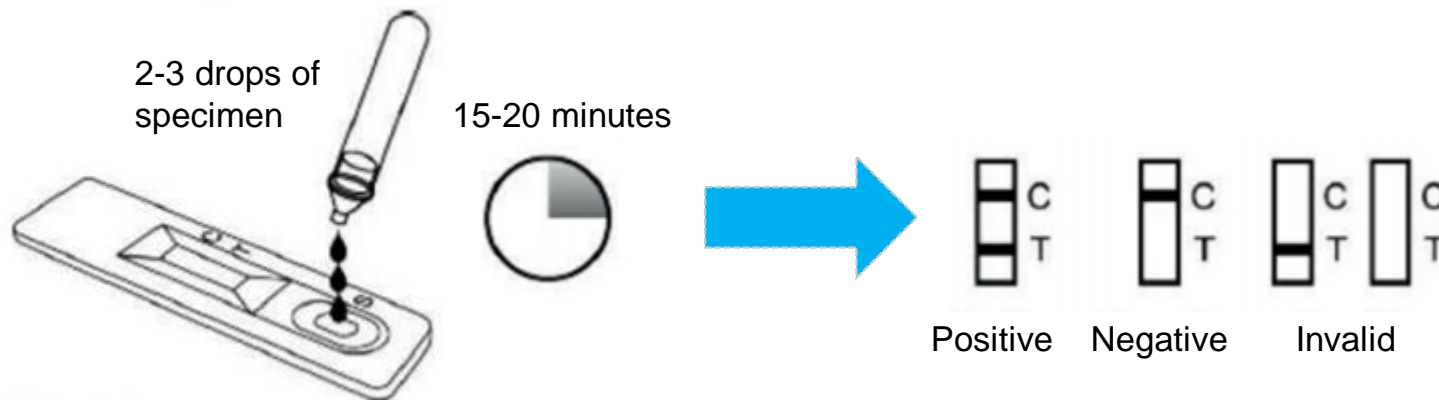

Ethically Priced To Be As Affordable As Possible

To order...

TEST KIT
EXPRESS



www.testkitexpress.co.uk



Call now on 01749 684085

Ethically Priced To Be As Affordable As Possible